


DECK MOUNT




CPAP Device

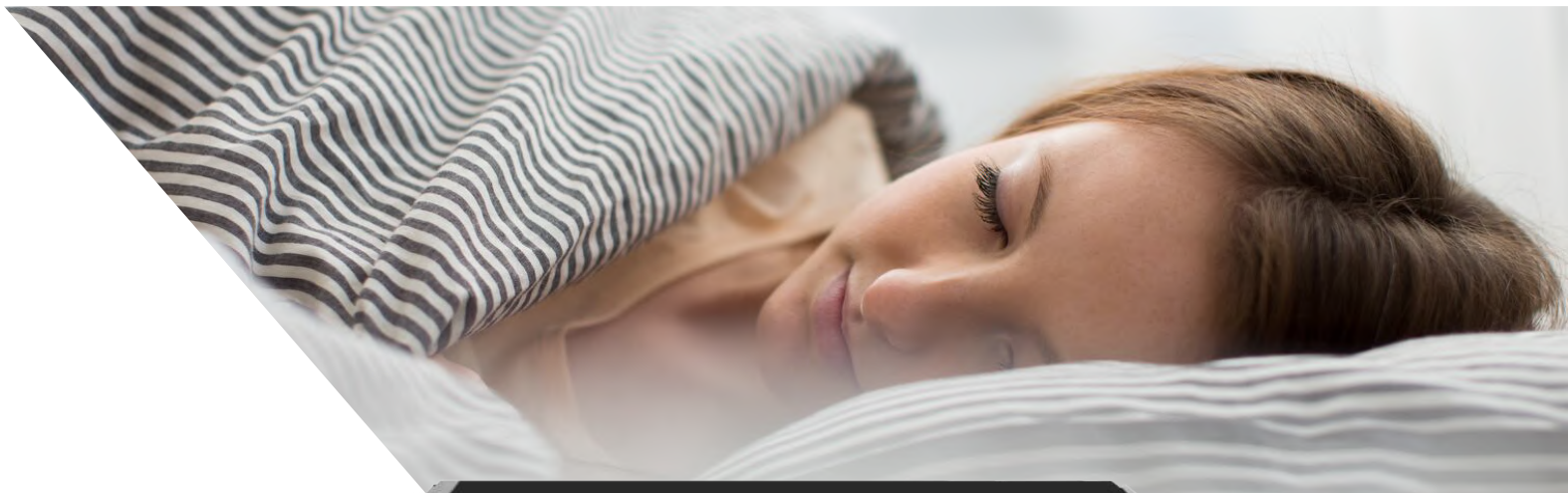

Harmony

Breathing easy, every night with CPAP

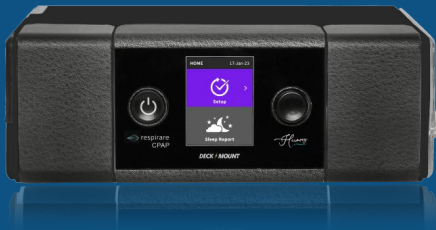

On Screen Report

Simplest way to get 365 days report in sleep report menu.



Gender Selection

Gender Selection helps to give a Soft/Standard response to the selected gender accordingly.

Technical Specification

Modes	CPAP, Auto CPAP
Operating/ Monitoring Range	4-20 CmH2O
Sound Pressure Level	<28 dBA with uncertainty of 2 dBA
Leak Compensation	Auto
Peak Flow	220 LPM
Ramp Time	Off, Auto, 5-45 mins
A-Flex (Expiratory Pressure relief)	On/ Off, 1-3 levels
Gender Selection	Male/ Female
Humidification Level	Off, 1-5 level (70°C max)
Display	2.8 in LCD
Data Management	1/7/31 days, 6 months, 1 year of detailed data.
Dimensions (LxWxH)	290x160x110 mm (Device with Humidifier)
Weight	1800 g (Device with Humidifier)
Power Adapter	Input Range : 100-240V, 50/60 Hz, 12V/ 60Watt, Over voltage protection Auto Controlled Temp variation
Battery Backup	External (2-8 Hrs)
Operating Temperature/ Humidity/ Altitude	+5 °c to +35 °c/ 10-95% relative humidity, non-condensing/sea level to 2,591 m, air pressure range 735 hpa to 1013 hpa
Storage and Transport Temperature/Humidity	-20 °c to +60 °c / 5-95% relative humidity, non-condensing
Supplement Oxygen	Maximum Flow: 4L/min (CPAP)
Standard Air Filter	Pollen Filter (Washable)
Tubing	Standard (19 mm)
Air Outlet	22 mm Air Outlet (STD)
Warranty	2 Years (Repair/ Replacement)*
Conforms to Standards	IEC 60601-1-1 IEC 60601-1-2 IEC 80601-2-12 CDSCO Approved
	Class 2 (double insulated, Type BF, Ingress Protection IP22)

*Term & Conditions Apply

Key Features

A-Flex

The A-Flex is a comfort feature that will lower the pressure when you exhale and slowly increases the pressure when you inhale and will be available in Auto CPAP mode.

Integrated Humidification

Heated Humidification makes therapy more comfortable and can reduce side effects like a dry nose and throat.

Mobility, Portability

Easy to carry the device from one location to another.


Product Code

Respirare Auto CPAP with Humidifier **VT 40 D**

Note: Technical specification may change without notice. Always refer to your device user manual.

Manufactured by:

Deck Mount Electronics Pvt. Ltd 260, Phase IV, Udyog Vihar, Sector 18, Gurugram, Haryana 122015

 1800 30 92 499

 Made in India

 wecare@deckmount.in

 www.deckmount.in